Dubliersilikon transparent

Duplicating silicone transparent

Produktinformation

Product information

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General product information

Addition-vulcanizing casting silicone for manufacturing negative molds in the production of earmolds and for the control molds.

Duplicating silicone, transparent is used primarily for processing UV-polymerizates or for producing control molds.

Final degree of hardness approx. 19–21 shore.


Delivery units Duplicating silicone, transparent

2 x 1 liter	REF 41000
2 x 5 liter	REF 41001
8 x double cartridges	REF 41003

Product information

Stability / Batch: See tub / cartridge

Storage: Protect against: frost, heat, UV-radiation / sunlight. Keep container tightly closed in a cool, well-ventilated place.

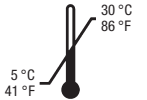Indications

Component A – catalyst

Component B – base

- Avoid any contact between the components A and B before mixing them.
- Close the receptacle tightly after manual use.
- Avoid direct skin contact with the liquid material! Before handling the material protect your hands by wearing suitable gloves.
- Cross-linked materials are chemically resistant – avoid stains on your clothes.

Notes for processing

Before processing the material it must have reached a temperature of $23\text{ °C} \pm 2\text{ °C}$ ($73.4\text{ °C} \pm 35.6\text{ °F}$) at a relative air humidity of $50\% \pm 5\%$.

Processing time: 4–6 minutes

Curing time without pressure pot: 15–30 minutes

Curing time with pressure pot: 10–15 minutes at 40 °C (104 °F) and 2–5 bar

Lower temperatures prolong and higher temperatures reduce time of processing and curing.

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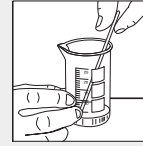
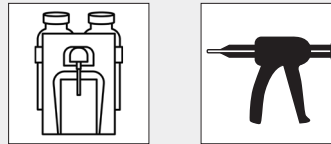
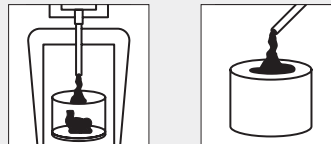
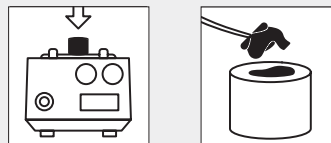
Processing

Manual mixing

Mix the two components A and B (mixing ratio 1:1) until material has homogeneous color. Approx. 40 – 50 g of duplicating silicone is needed for one negative mold.

Automatic mixing

The material can be processed with the silicone-mixing device SIDOMIX II (REF 54561) or with an inject pistol (REF 25100). Mix the components A and B bubble-free in the mixing canula (mixing ratio 1:1). Do not use mixing canula (REF 25151 ff or REF 61402 ff) if it contains hardened material.

Filling

When filling the material into the investment mold make sure that the tip of the mixing canula always stays in the material to avoid the formation of bubbles.

Polymerization

Afterwards let the material cure at room temperature or put the filled investment mold in the pressure pot (REF 54200/49800) for 10 – 15 minutes at 2 – 5 bar for a bubble-free polymerization. After the polymerization of the duplicating silicone, use the impression removing instrument (REF 46700) to carefully demold the impression. Before continuing with processing the impression, check the negative mold for wax remainders and remove them.

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7908609_2022_05

Please follow the product information and the material safety data sheet.